



WHAT A CATCH!

QuEChERS Cartridges for Pesticide Residue Analysis

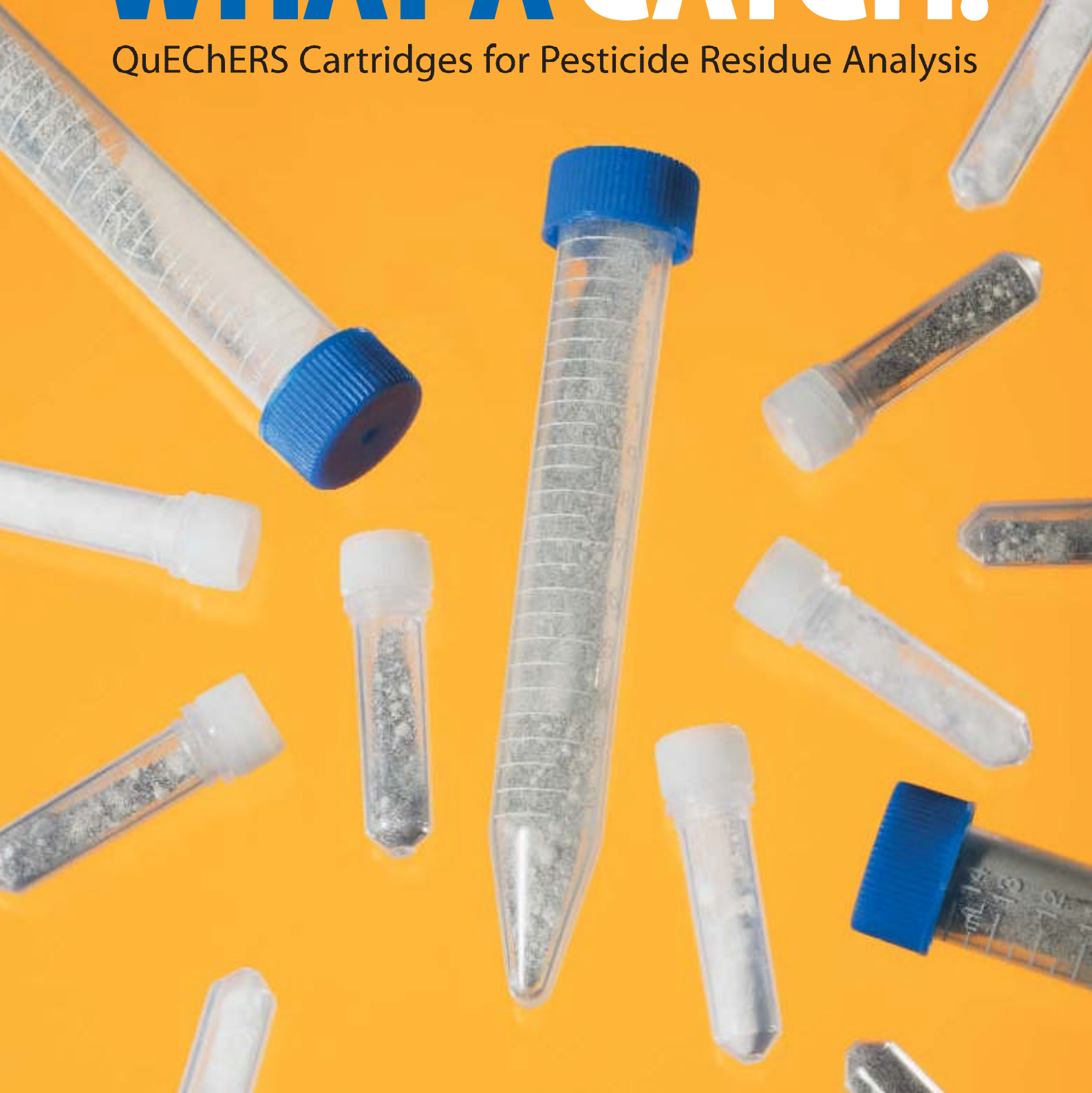

Chromatography Products

www.restek.com

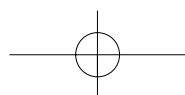

QuEChERS SPE Cartridges

Quick, Easy, Cheap, Effective, Rugged, and Safe, the QuEChERS (“catchers”) method developed by the USDA Eastern Regional Research Center, has become very popular for extraction and clean-up of pesticide residue samples.

- Fast, simple sample cleanup
- Fourfold increases in sample throughput
- Fourfold decreases in material cost
- Variety of formats, to meet all needs
- SPE dispersive technique

Multiple sorbents:

- PSA* - removes sugar and fatty acids
- MgSO₄ - removes water
- C18 - removes nonpolar interferences
- Carbon - removes pigments & sterols
- Custom combinations for your specific needs



Custom products are available by quote request.

Contact your local distributor in your country or visit www.restek.com/quechers



Lit. Cat.# 580168-INT

© 2007 Restek Corporation.



Free Sample Packs Available!

To receive your free 3-pack or 5-pack, add -248 to the item number.
(One sample per customer.)

QuEChers SPE Cartridges	qty.	cat#
QuEChERS SPE 2mL Micro-Centrifuge Cartridge Packed with 150mg Magnesium Sulfate and 50mg PSA	100-pk.	26124
QuEChERS SPE 2mL Micro-Centrifuge Cartridge Packed with 150mg Magnesium Sulfate, 50mg PSA, and 50mg Graphitized Carbon	100-pk.	26123
QuEChERS SPE 2mL Micro-Centrifuge Cartridge Packed with 150mg Magnesium Sulfate, 50mg PSA, and 50mg C18	100-pk.	26125
QuEChERS SPE 15mL Centrifuge Cartridge Packed with 900mg Magnesium Sulfate, 300mg PSA, and 150mg Graphitized Carbon	50-pk.	26126
QuEChERS SPE 6mL SPE Cartridge Packed with 200mg Graphitized Carbon and 400mg PSA, Teflon® Frits	30-pk.	26127
QuEChERS SPE 6mL SPE Cartridge Packed with 250mg Graphitized Carbon and 500mg PSA, Teflon® Frits	30-pk.	26128
QuEChERS SPE 6mL SPE Cartridge Packed with 500mg Graphitized Carbon and 500mg PSA, Teflon® Frits	30-pk.	26129

*PSA: Primary and secondary amine exchange material.

Available from:

LECO Australia Pty Ltd
PO Box 6006
Baulkham Hills BC NSW 2153

Ph: 61 2 9894 5955
Fx: 61 2 9894 5247

Email: australia@leco.com

